


# Copper Pigtails

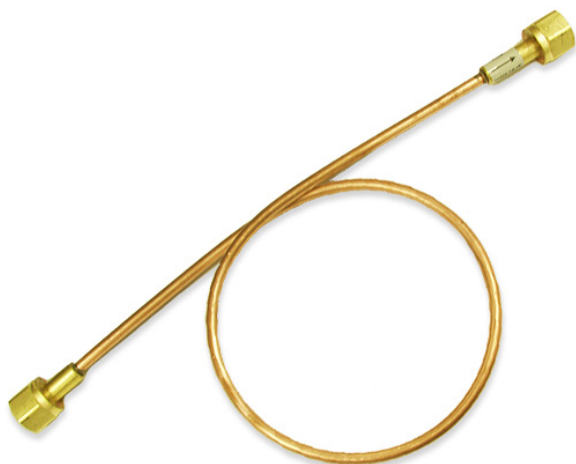

## Features:

- Burst pressure rating at room temperature: 12,000 psig
- Maximum 3000 psig working pressure on tubing
- All pigtails 100% tested before shipping and labelled for gas service and pressure rating
- Standard length 36" (3')
- Includes a sticker that indicates the direction of the flow

## General Specifications:

The Amico copper pigtails shall have a 5/16" x 0.065" (8 mm x 1.7 mm) wall thickness, copper tube. The Amico copper pigtails shall come with a check valve equipped to one of the ends.


Note: Burst pressure ratings for pigtails are determined at room temperature with the hose in a straight line. A safety factor of 4:1 should be used for normal applications. Impulse and shock pressure applications may require a higher safety factor.

Note: Do not apply heat to pigtails.

All cylinder connections machined to CGA specifications for applicable gases.

Amico products comply with NFPA 99 and CSA Z7396.1.

# Technical Specifications

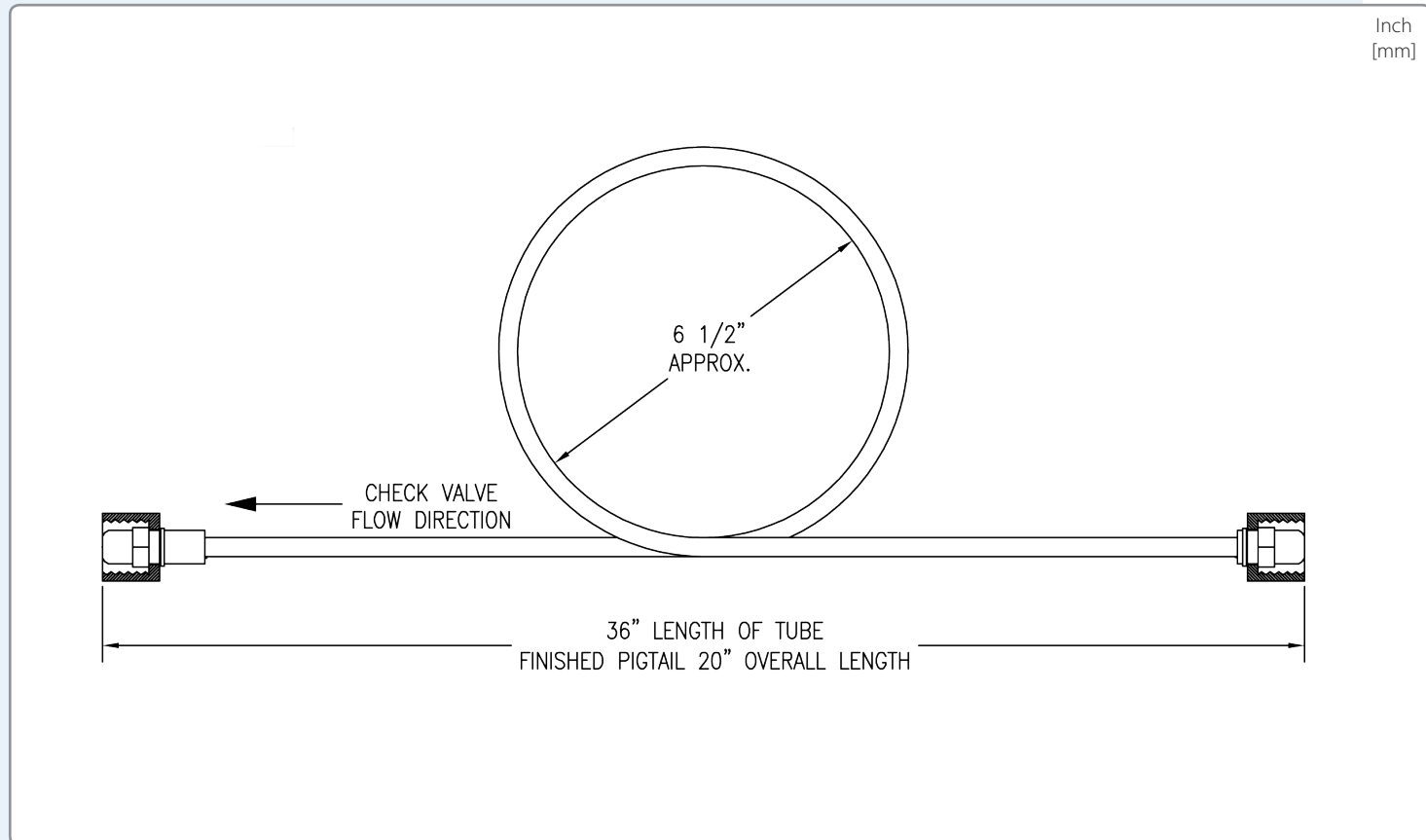

## Model Numbers

### **M-X-HB-PTC-XXX**

**The XXX defines the Gas:**

Nitrogen	= NIT
Medical Air	= AIR
Carbon Dioxide	= CO2
Nitrous Oxide	= N2O
Oxygen	= OXY